


Single-cell RNA-seq Preparation and Best Practices



General considerations to increase the likelihood of high quality data

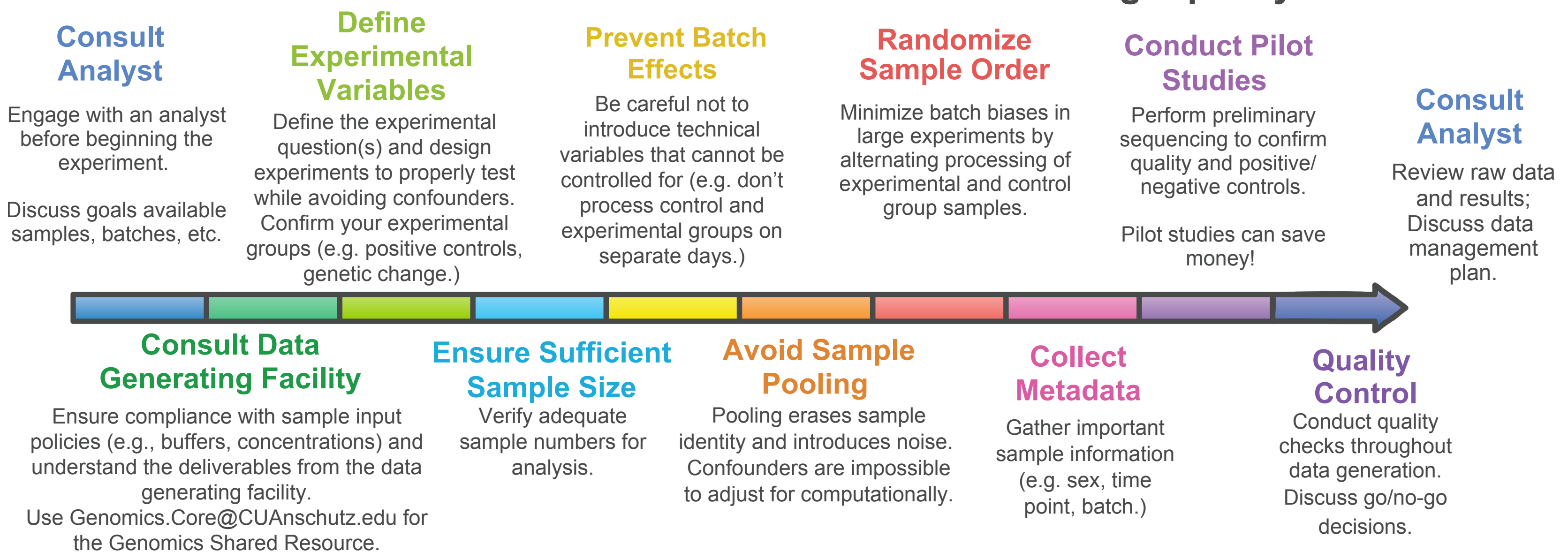

scRNA-seq captures and sequences the transcriptome of individual cells

Used to explore transcriptional changes between and within cell populations, either at baseline or after treatment, gene KO, etc.

Define the objective

Outline the biological question, cells of interest, and the primary readouts before preparing samples.


Plan cell numbers & enrich strategically

Adjust the input cell number to the rarest population desired. Use positive or negative selection (e.g. FACS/microfluidics, bead- or column-based kits) to enrich target populations.

Consider reducing batch effects (and cost) with sample multiplexing

Current options include on-chip multiplexing (OCM - recommended), fixed RNA-profiling (Flex), hashtag oligos (HTO), cell multiplexing oligos (CMO), and genotype-based methods.

Confirm the approach is compatible with the cell barcode chemistry, and understand that some methods can result in experimental failures.

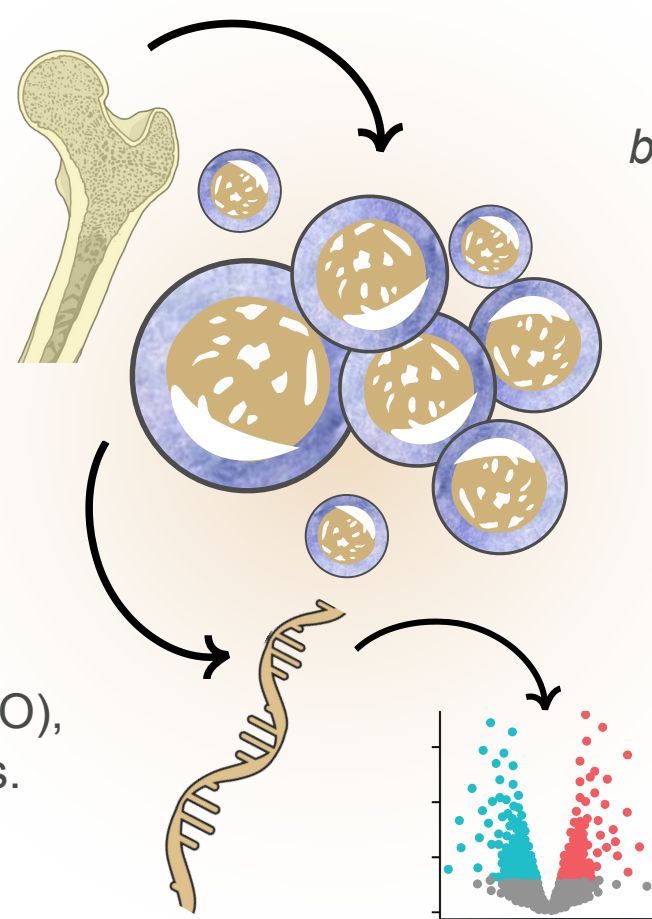

Sample viability drives data quality

Aim for >80-90% viability per sample and optimize each step *before* sequencing. Cell membrane integrity is critical to produce robust and interpretable results.

Deplete dead/dying cells and remove debris

Use positive or negative selection (e.g. FACS/microfluidics, EasySep/Annexin V) to remove dead cells.

Clear debris and cell aggregates with extra large volume washes and cell strainer/filters (e.g. MACS Smart Strainer.)

Calculate cell concentrations accurately

Avoid skewed target recovery by calculating stock concentrations with the *total* cell count (not just viable cells.)

Count the final suspension 2-3x; individual counts can vary as much as 25%.

Before sequencing

- ☐ Discussed experiment with an analyst & sequencing facility
- ☐ Defined target cells & questions
- ☐ Finalized sample multiplexing strategy
- ☐ Budgeted for throughput and sequencing depth
- ☐ Tissue dissociation is optimized
- ☐ Samples measure >80-90% viable
- ☐ Depleted dead / dying cells
- ☐ Removed debris / cell aggregates
- ☐ Calculated 2-3 cell counts that are <25% variable

After sequencing

- ☐ Remove empty droplets / barcodes
- ☐ Identify doublets / multiplets
- ☐ Filter low quality cells and samples
- ☐ Measure ambient RNA (correct if needed)

Some analyses will require experience with a coding environment (R/Python)

Library Type	3'/5' scRNA-seq	Fixed (Flex) scRNA-seq	Long read scRNA-seq
Protocol	Droplet-based UMIs (10x Genomics)	Droplet-based UMIs (10x Genomics)	Plate-based full-length transcript (SmartSeq)
Gene Features	Whole-transcriptome (3' or 5' mRNA)	Probe-based (option to include custom probes)	Whole-transcriptome (includes non-coding molecules)
Cells Captured	<20k cells / sample	<20k cells / sample	100-1k cells / sample
Sequencing Depth	50-100k reads / cell	10-40k reads / cell	>1M reads / cell
More reads may be needed for lowly expressed genes or RNA-rich samples (e.g. cell lines)			
Multimodal Options	ATAC, CITE/protein, TCR/BCR, CRISPR	CITE/protein, CRISPR	Generally not compatible
Common Analyses (beyond clustering)	Differential expression Differential abundance Trajectory inference RNA velocity VDJ clonotypes CNV/CNA inference Cell-cell communication	Differential expression Differential abundance Trajectory inference CNV/CNA inference Cell-cell communication	Alternative splicing & isoform detection Differential expression Differential abundance Trajectory inference CNV/CNA inference SNP/variant calling Cell-cell communication
≥2 biological replicates are recommended (more for <i>in vivo</i>)			

Contact the Biostatistics & Bioinformatics Shared Resource (BBSR) for cancer-related analysis

James.Costello@CUAnschutz.edu, Bioinformatics@CUAnschutz.edu



Contact CUBE for all non-cancer-related projects



scRNA-seq v2.2